Feline Leukemia Virus Antigen Test Kit

For veterinary use only.

English version

SNAP™ FeLV

SNAP FeLV Antigen test is an enzyme-linked immunosorbent assay (ELISA) for the detection of feline leukemia virus (FeLV) antigen in feline serum, plasma or whole blood. The detection of the FeLV group-specific viral antigen (p27) is diagnostic for FeLV infection. The SNAP FeLV Antigen assay utilizes monoclonal antibodies to p27 and positive controls.

The conjugate contains enzyme-conjugated antibody to p27. Upon mixing the conjugate and the test sample, conjugated monoclonal antibody will bind p27 antigen (if present). The sample/conjugate mixture is then added to the SNAP device and flows across the spotted matrix. The matrix-bound p27 antibody (FeLV spot) will capture the p27-conjugated antibody complex. The device is then activated, releasing wash and substrate reagents stored within the device.

Color development in the FeLV Ag sample spot indicates the presence of FeLV antigen.

Precautions and warnings

- All wastes should be properly decontaminated prior to disposal.
- Do not mix components from kits with different serial numbers.
- Do not use a SNAP device that has been activated prior to the addition of sample.
- Refer to country specific Safety Data Sheet for regional hazard identification.

Storage

Store devices and test reagents at 2–8°C until the expiration date. All components must be at room temperature (18–25°C) before running the test—**do not heat**.

Kit components

Item	Reagents	Quantity
1	1 bottle anti-FeLV: HRPO conjugate (Preserved with ProClin™ 150)	8.0 mL
2	SNAP device	15 or 30
Reagents contained in each device:		
	Wash solution (preserved with ProClin™ 150)	0.4 mL
	Substrate solution	0.6 mL
Other Components: transfer pipettes, sample tubes, and reagent rack		

WARNING

Conjugate - H317 / H402 / H412 / P261 / P272 / P273 / P280 / P302+P350 / P363 / P333+P313 / P501: May cause an allergic skin reaction. Harmful to aquatic life. Harmful to aquatic life with long lasting effects. Avoid breathing mist/vapors. Contaminated work clothing must not be allowed out of the workplace. Avoid release to the environment. Wear protective gloves. IF ON SKIN: Wash with plenty of water. Wash contaminated clothing before reuse. If skin irritation or rash occurs: Get medical advice/attention. Dispose of contents/container in accordance with local/regional/national/international regulations.

Wash Solution - H317 / H319 / H402 / H412 / P261 / P264 / P272 / P273 / P280 / P305+P351+P338 / P337+P313 / P333+P313 / P363 / P501: May cause an allergic skin reaction. Causes serious eye irritation. Harmful to aquatic life. Harmful to aquatic life with long lasting effects. Avoid breathing mist/vapours. Wash thoroughly after handling. Contaminated work clothing must not be allowed out of the workplace. Avoid release to the environment. Wear eye protection/face protection. Wear protective gloves. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. If eye irritation persists: Get medical advice/attention. If skin irritation or rash occurs: Get medical advice/attention. Take off contaminated clothing and wash it before reuse. Dispose of contents/container in accordance with local/regional/national/international regulations.

Sample information

- Samples must be at room temperature (18–25°C) before beginning the test procedure.
- Serum, plasma or anti-coagulated whole blood (e.g., EDTA, heparin), either fresh or stored at 2–8°C for up to one week, can be used.
- For longer storage, serum or plasma can be frozen (-20°C or colder) and then recentrifuged before use.
- Hemolyzed or lipemic samples will not affect results.

Test procedure

1. Allow all components to equilibrate at room temperature (18–25°C) for 30 minutes before use. **Do not heat**.

2. Using the pipette provided, dispense **3 drops of sample** into a new sample tube.

3. Holding the bottle vertical, add **4 drops of conjugate** to the sample tube.


4. Cap the sample tube and mix it thoroughly by **inverting it 3–5 times**.

5. Place the device on a horizontal surface. Add the entire contents of the sample tube to the sample well, being careful not to splash the contents outside of the sample well.

 The sample will flow across the result window, reaching the activation circle in 30–60 seconds. Some sample may remain in the sample well.

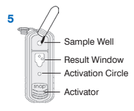6. When color FIRST appears in the activation circle, push the activator firmly until it is flush with the device body.

 Note: Some samples may not flow to the activation circle within 60 seconds and, therefore, the circle may not turn color. In this case, press the activator after the sample has flowed across the result window.


7. Test results must be read at **10 minutes** from the time of activation.

 Note: The positive control may develop sooner, but results are not complete until 10 minutes.

Interpreting test results

Positive result

Any color development in the sample spots indicates the presence of FeLV antigen in the sample.

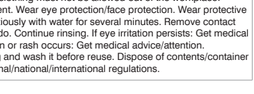Negative result

Only the positive control spot develops color.

Invalid results

- **Background**—If the sample is allowed to flow past the activation circle, background color may result. Some background color is normal. However, if colored background obscures the test result, repeat the test.
- **No color development**—If the positive control does not develop color, repeat the test.

Symbol Descriptions

	Use by	EC REP	Authorized representative in the European Community
LOT	Batch code (Lot)		Consult instructions for use
SN	Serial number	IVD	In vitro diagnostic
	Temperature limitation		Date of manufacture
	Manufacturer		
REF	Catalog number		

IDEXX Technical Support
USA/Canada: 1 800 248 2483 • idexx.com
Australia: 1300 44 33 99 • idexx.com.au
Europe: idexx.eu

VLN/PCN: 313/5028.04

SNAP is a trademark or registered trademark of IDEXX Laboratories, Inc. or its affiliates in the United States and/or other countries.

Patent information: idexx.com/patents

© 2025 IDEXX Laboratories, Inc. All rights reserved. • 06-0021123-00

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© 2025 IDEXX Laboratories, Inc. Všetky práva vyhradené.
IDEXX B.V.